



-80°C
-40°C
-20°C
RT


SAB
Series

BiO.Li.X!


SAB-SERIES - AUTOMATED BIO-LIBRARIES

The SAB-Series combines sample safety, sample density, and power efficiency design principles of the BiOLiX! products, but also accommodates the very specific operational needs of collections with diverse and unpredictable tube and rack types. The SAB also provides an extremely easy pathway to transition manual freezer repositories into automation operations.

Biobanking is becoming more and more complex with growing sample types, methods and organizations. Centralization of collections and services is often a key component of biobanking operations. Large sample collections with diverse sample types and non-standardized labware present a challenge to operational and cost efficiency. LiCONiC Instruments is a leader in dedicated, application driven sample storage automation. The SAB series of the BiOLiX! Bio-Library Sample Storage product family specifically addresses the need of this rapidly growing sample management sector.

Energy Efficiency

The SAB-Series delivers an industry leading return on investment by all standards. The combination of a highly space optimized storage layout, unique chest freezer design and our proprietary refrigeration design, reduces power consumption to a fraction of what is known in the market today.

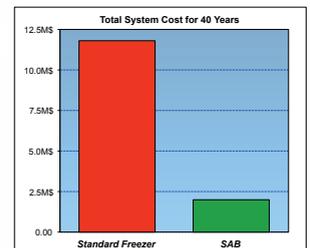

MAIN FEATURES

- Temperature range -25°C to -80°C
- Semi-automated / high throughput
- Large sample capacity
- Exceptional frost control
- Individual sample tracking and picking option

SAB-SERIES PRODUCT OVERVIEW

The SAB-ULT-Series is designed for biobanking applications where the primary objective of the storage is long term sample safety and efficiency for storage processes and capacity. This series offers the same advantages as other LiCONiC store lines but caters to the needs of large legacy collections and long term collection arrangements, where odd and variable sized containers and rack sizes are frequently encountered. The SAB-ULT-Series simplifies the sample storage/retrieval process and ensures that sample integrity is maintained in perpetuity.

The SAB storage area is flexibly configurable to accommodate many container varieties simultaneously and to be changeable over time to adapt to new container types and proportions in storage.

Transitioning sample collections from legacy storage arrangements to the automated environment of the SAB is easy with fully configurable storage cassette sizes.


Ideal for legacy cryobox collections -
No reformatting required.



The SAB-ULT-Series is designed for the unique challenges of large biobank sample collections.

Examples of appropriate applications:

- Genetic or disease based population cohort studies
- Forensic samples studies
- Cellomic samples
- Genomic samples
- Transitioning large collections from manual to automated storage
- Core facility receiving many odd or unpredictable sample storage
- Container varieties

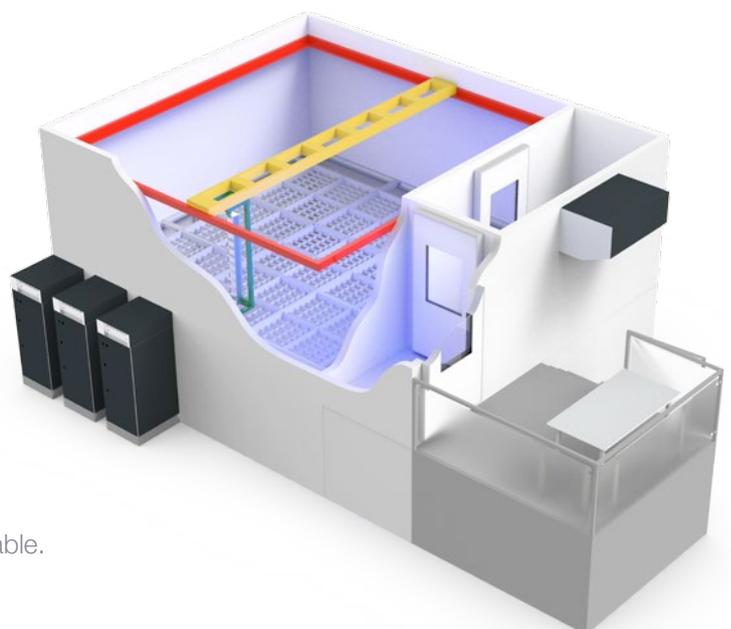

Full rack or picked, individual sample retrieval options available.

SAB CAPACITIES

BiOLiX! SAB-Series examples			SAB46k0-ULT		SAB98k0-ULT		SAB220k5-ULT	
Labware	Plate height	Pitch	Levels	Samples	Levels	Samples	Levels	Samples
FluidX 0.3 ml Screw Cap	19.1 mm	25 mm	50	4'147'200	78	8'626'176	78	19'408'896
Micronic 0.5 ml Screw Cap	28.5 mm	35 mm	36	2'985'984	56	6'193'152	56	13'934'592
LVL LX1000	45.6 mm	52 mm	24	1'990'656	37	4'091'904	37	9'206'784
Thermo Matrix 1.0ml Screw Cap	54.6 mm	61 mm	20	1'658'880	32	3'538'944	32	7'962'624
Greiner Bio One 2 ml	48.0 mm	54 mm	23	953'856	36	1'990'656	36	4'478'976




SAB STORAGE CONDITIONS

The SAB-Series features climate control design found across the BiOLiX! product family. These features result in high density, temperature stable, and power efficient sample storage. As with the SAB-Series, single, chest style storage combines with LiCONiC's proprietary dual temperature

zone design to provide temperature stability, excellent humidity control and easy service access to robotic components.


LiCONiC's unique chest freezer design naturally creates numerous benefits. Two of the most important are:

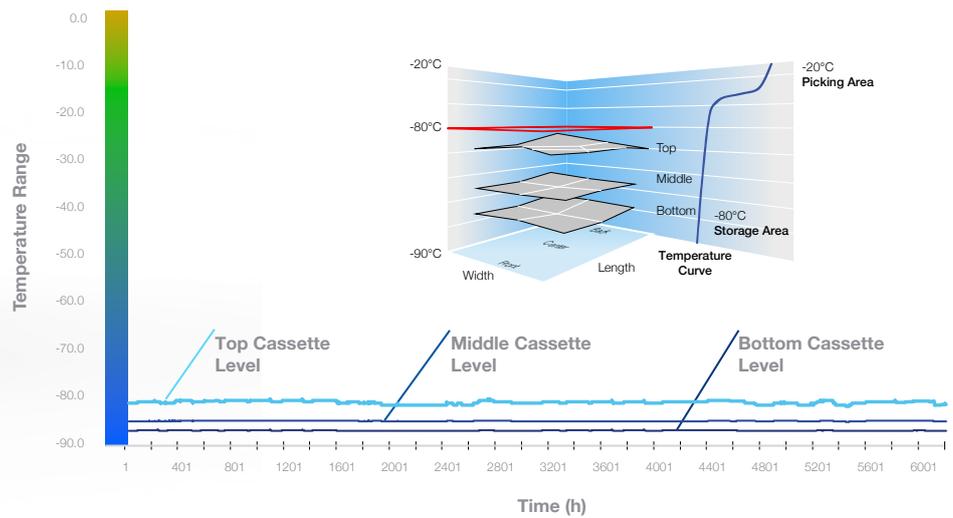

Two climate zones, one for active components at -20°C and another for passive elements at -80°C. Proven, stable horizontal temperature layers within the ULT storage area.

SPACE EFFICIENCY

The SAB-Series Automated BioLibraries makes best use of given room constraints by its "built to space" modular design. System dimensions are optimised to the individual application's needs and workflow requirements.

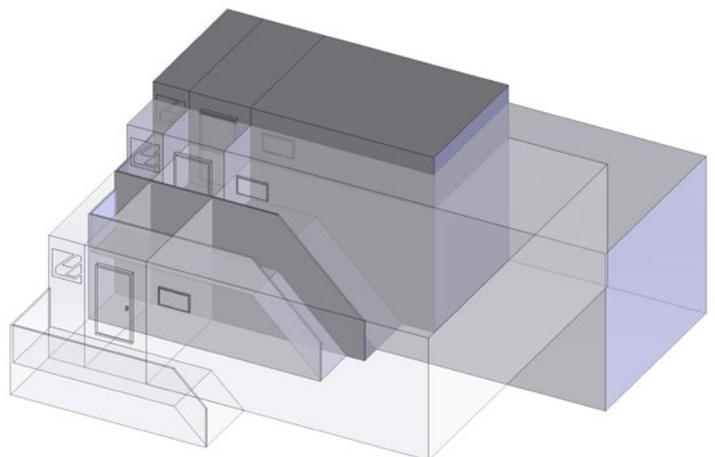




INTUITIVE SOFTWARE

- Graphical Job Creation for Import/Export of Samples
- Easy Search of Samples by Barcode, Name, Data etc.
- Complete Audit Trail
- Automatic Inventory Defragmentation
- Automatic Data Backup
- Integration to LIMS Database
- Remote Support
- Mobile Device Control
- User Access Security

OPTIONS & ACCESSORIES

- Barcode Reader
- Flexible Labware Support
- Configurable to Room Constraints
- Tablet Interface
- External Database Integration
- Remote Monitoring
- Internal Vision System
- Multi - Tube / Vial Picker
- LN2 - Backup Option
- 2D - BCR Rack Scanner
- Redundant Refrigeration
- Dynamically Configurable Storage Cassettes

SAB DIMENSIONS

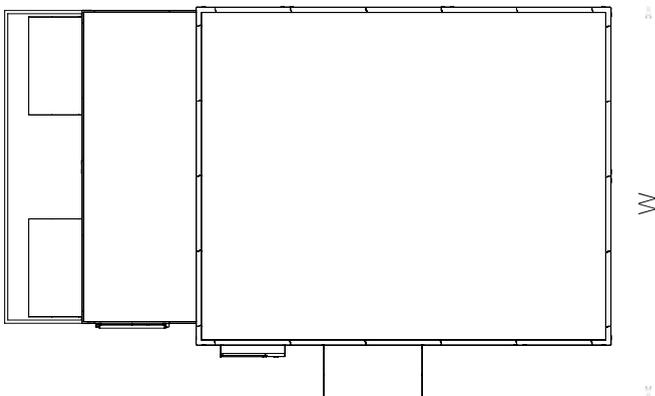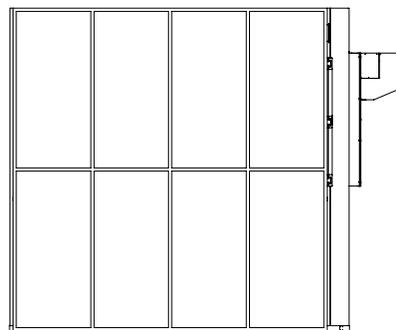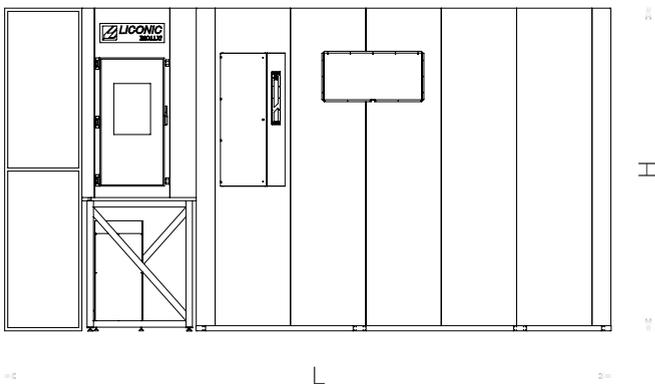

BiO.Li.X!		SAB Dimensions examples		
Type	L= Length	W= Width	H= Height	
SAB46k0-ULT	8045 mm	6092 mm	3530 mm	
SAB98k0-ULT	9410 mm	6092 mm	5030 mm	
SAB220k5-ULT	15845 mm	6092 mm	5030 mm	

WORLDWIDE PRESENCE



CORPORATE **LiCONiC AG**

LiCONiC AG
Industriestrasse 8-12
9493 Mauren
Principality of Liechtenstein

Phone: +423 340 5000
Fax: +423 37 35 35 9
E-Mail: info@liconic.com
sales@liconic.com
service@liconic.com

LiCONiC Services **Deutschland GmbH**

Bahnhofsplatz 3
D-56410 Montabaur
Germany

Phone: +49 2602 997 4557
Fax: +49 2602 997 4559
E-Mail: info.de@liconic.com
sales.de@liconic.com
service.de@liconic.com

LiCONiC UK Ltd

BioHub
Alderley Park
Macclesfield
Cheshire, SK10 4TG
United Kingdom

Phone: +44 (0) 1618 752 542
+44 (0) 1618 752 543
Fax: +44 (0) 1618 752 540
E-Mail: info.uk@liconic.com
sales.uk@liconic.com
service.uk@liconic.com

LiCONiC US, Inc.

21-F Olympia Ave
Woburn, MA 01801
USA

Phone: +1 (781) 309-3340
Fax: +1 (781) 537-6173
E-Mail: info.usa@liconic.com
sales.usa@liconic.com
service.usa@liconic.com

LiCONiC Japan K.K.

2-2-15 Daiya Bldg 2F
Hamamatsu-cho, Minato-ku
Tokyo Japan 105-0013

Phone: +81 3 6841 0707
E-Mail: info.japan@liconic.com
sales.japan@liconic.com
service.japan@liconic.com

LiCONiC China Co. Ltd

1-B102, No.3891, Jinxiu Rd
Shanghai (201204),
R.P.China

Phone: +86 186-2156-0603
E-Mail: info.china@liconic.com
sales.china@liconic.com
service.china@liconic.com

Asia:

Phone: +1 (760) 522-2455
E-Mail: info.asia@liconic.com
sales.asia@liconic.com
service.asia@liconic.com



Information within this document is subject to be changed without prior notice. For details please contact our applications at info@liconic.com. 20190411



www.liconic.com